

ORIPLAST®



Secretion suction

GENTLY CLEAR

Secretion suction

AIRWAY CLEANING

Diagnostic and therapeutic secretion suction

Products for secretion suction are used both diagnostically and therapeutically: On the one hand, they enable the hygienic collection of samples from the ear, nose and throat area, on the other hand, the rapid removal of liquids that expose the airways. The latter plays a role, especially in newborns, in order to enable independent and easier breathing.


The high-quality transparent collection container has an optimal capacity for the manually drained secretion.

Collection container

The long, double-barreled tubing makes it possible to identify a stop or flow of secretion thanks to its transparent properties.

Tubing



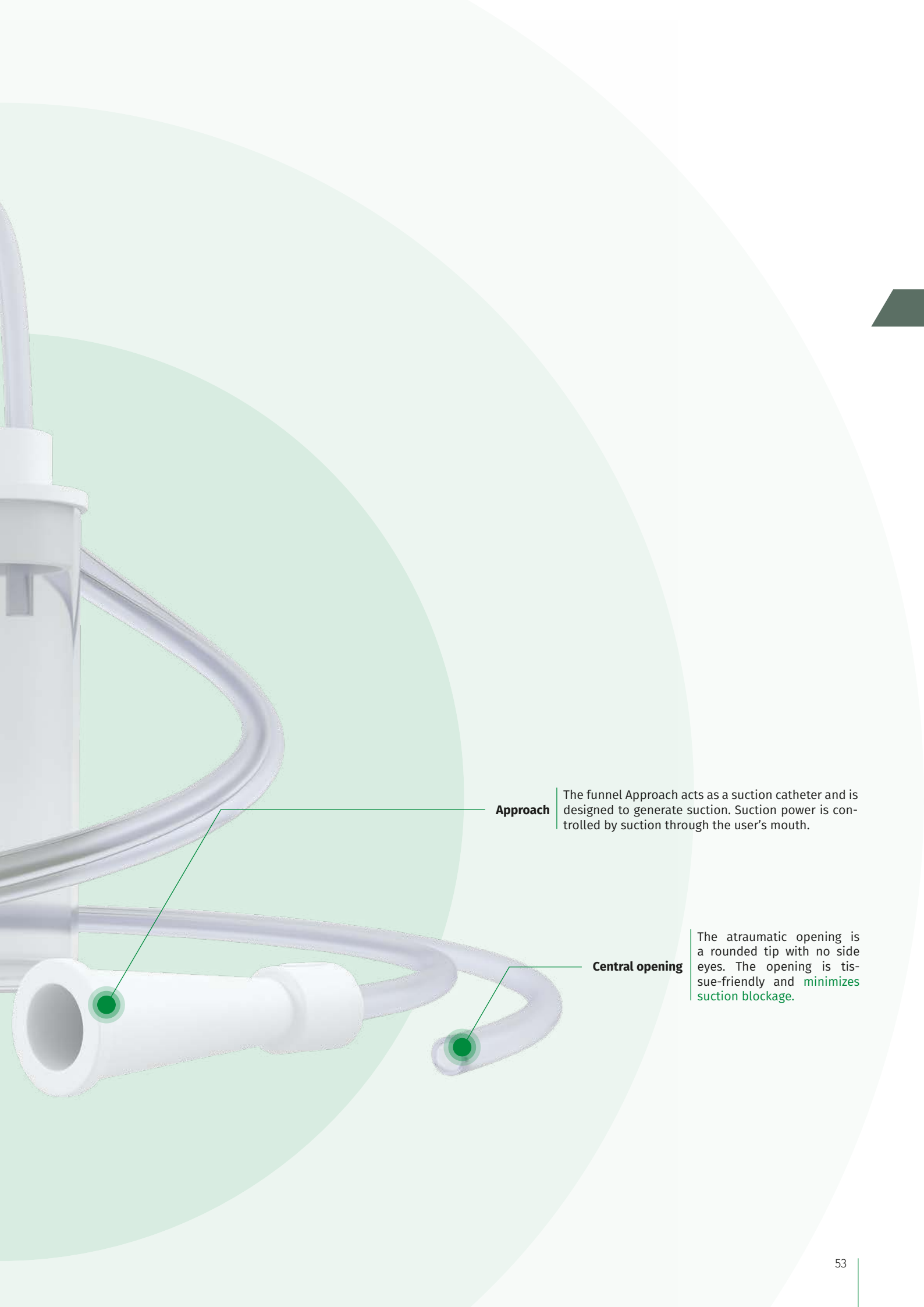

Approach

The funnel Approach acts as a suction catheter and is designed to generate suction. Suction power is controlled by suction through the user's mouth.

Central opening

The atraumatic opening is a rounded tip with no side eyes. The opening is tissue-friendly and minimizes suction blockage.

Suction device

The latex-free mucus aspirator is specially designed for use with newborns. The suction power is controlled by suction through the helper's mouth. The Tracheal Suction Set is used to collect tracheal secretion for laboratory testing. After collection, the lid with the suction Approach s is removed and replaced with the sterile stopper provided.



Slime extractor

- For newborns
- Volume 25ml
- Atraumatic tubing end

eClass 34-24-13-90



STERILE EO

Item No.	PZN	Tubing length (mm)	Approach	FG	DEHP	PVC	PU
VK00941	17488738	2 x 450	Funnel	12			50
VK00942	17488744	2 x 450	Funnel	8			50



Tracheal suction set

- for collecting tracheal secretions
- Volume 13.5ml
- CH 14 funnel Approach
- enclosed sealing plug

eClass 34-24-13-90



STERILE EO

Item No.	PZN	Tubing length (mm)	Connector	DEHP	PVC	PU
VK00944	17488750	80 + 100	Step connector			50
VK00945	17488767	2 x 80	Fingertip			50